



Hillrom™



# Hillrom™ P290 Surface

PROVIDING A SAFE  
SKIN ENVIRONMENT  
FOR YOUR PATIENTS


Pressure ulcers are costly and prevalent. That's why the P290 mattress overlay or mattress replacement system from Hillrom was designed to help prevent and treat pressure ulcers in low to high risk adult patients. The P290 surface provides an adjustable level of immersion and envelopment in both Alternating and Continuous Low Pressure therapy modes.

This therapeutic approach reduces compression of vulnerable tissues and is aligned with The NPIAP/EPUAP/PPPIA International Guidelines which advise clinicians to change from a support surface to a specialty support surface (i.e. alternating pressure air mattress, low-air-loss mattress or air fluidized bed) when the individual<sup>1</sup>:

- Cannot be positioned off the existing pressure ulcer
- Has a pressure ulcer on two or more turning surfaces which limits repositioning positions
- Has a pressure ulcer that fails to heal or that deteriorates despite appropriate comprehensive care
- Is at high risk for additional pressure ulcers
- Has undergone flap or graft surgery
- Is uncomfortable
- 'Bottoms out' (i.e. touches the bed frame) on the current support surface



## Hillrom™ P290 Surface OVERLAY AND MATTRESS REPLACEMENT THERAPY SYSTEM

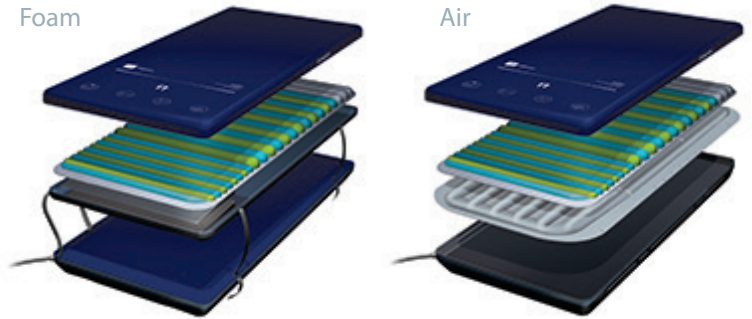

### FEATURES AND BENEFITS

Specifically designed for enhanced prevention and treatment of patients at risk or suffering from pressure ulcers.

A choice of alternating or continuous low pressure therapies to suit patient needs and clinical preference.

Adjustable cycle times and firmness setting for greater control over patient therapy and comfort.

Patient and caregiver safety features including automatic control panel lock-out, fault indicators and 2-hour transport mode.

Three pressure-redistributing zones provide additional protection for vulnerable body areas, including 5° sloped alternating heel zone.

P-MAX and PLUS inflation for additional support and safety during patient positioning and transfer activities.

Mattress replacement is available with either air or foam patient support layer.

X-ray sleeve to aid in bed x-rays and patient comfort.\*

Evacuation handles to assist with moving patients during an emergency.\*

\* Optional

### TECHNICAL SPECIFICATIONS

#### MATTRESS OVERLAY

Length	200 cm
Width	90 cm
Height	10 cm
Weight	6.1 kg

#### MATTRESS REPLACEMENT - FOAM

Length	200 cm
Width	90 cm
Height	17 cm
Weight	12 kg

#### MATTRESS REPLACEMENT - AIR

Length	200 cm
Width	90 cm
Height	17 cm
Weight	7.5 kg

#### CONTROL UNIT

Length	90 cm
Width	20 cm
Height	12 cm
Weight	2.2 kg

Safe Working Load	250 kg
-------------------	--------

Recommended Therapeutic Weight	30 to 140 kg
--------------------------------	--------------

#### Cover Material

Polyurethane coating on polyester material, 360° zipper, welded seams, low friction, two-way stretch, moisture vapour permeable, bacteriostatic, fungistatic and antimicrobial. Can be wiped and washed. Chlorine-based cleaning solutions may be used (≤1,000ppm).

**hillrom.com**

<sup>1</sup>EPUAP/NPIAP/PPPIA International Guidelines for the Prevention & Treatment of Pressure Ulcers, 2019.

The P290 Mattress helps to prevent and treat pressure ulcers in low- to high-risk adult patients, within the recommended patient weight limits of 30 to 140 kg, in order to achieve validated clinical performance in all the typical head of bed positions.

Class IIa

Notified Body: GMED CE0459

Manufacturer: Hill-Rom S.A.S. – BP 14 ZI du Talhouët 56330 Pluvigner France

This medical device is a regulated health product which, pursuant to such regulation bears a CE mark. Hill-Rom recommends that you carefully read the detailed instructions for safe and proper use included in the documents accompanying the medical devices. The personnel of healthcare establishments are responsible for the proper use and maintenance of these medical devices. Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.